

Patient:

DOB:

Age:

Sex:

Patient Identifiers:

Visit Number (FIN):

Client:

Physician:

ARUP Test Code: 3002528

Collection Date: 12/30/2021

Received in lab: 12/31/2021

Completion Date: 01/17/2022

COMMENT:

Clinical Information: 3 isolates of *Staphylococcus aureus* compared.

1,12115557

2,12117107

3,12116600

GROUPED INDISTINGUISHABLE ISOLATES:

Group A: 12115557, 12116600 and 12117107

GROUPS OR STRAINS NOT RELATED TO OTHER ISOLATES:

Group-A




Patient:

ARUP Accession: 21-364-101995

Bacterial Strain Typing by Next Generation Sequencing

Patient: | Date of Birth: | Sex: | Physician:
 Patient Identifiers: | Visit Number (FIN):

Dendrogram and Identity Matrix

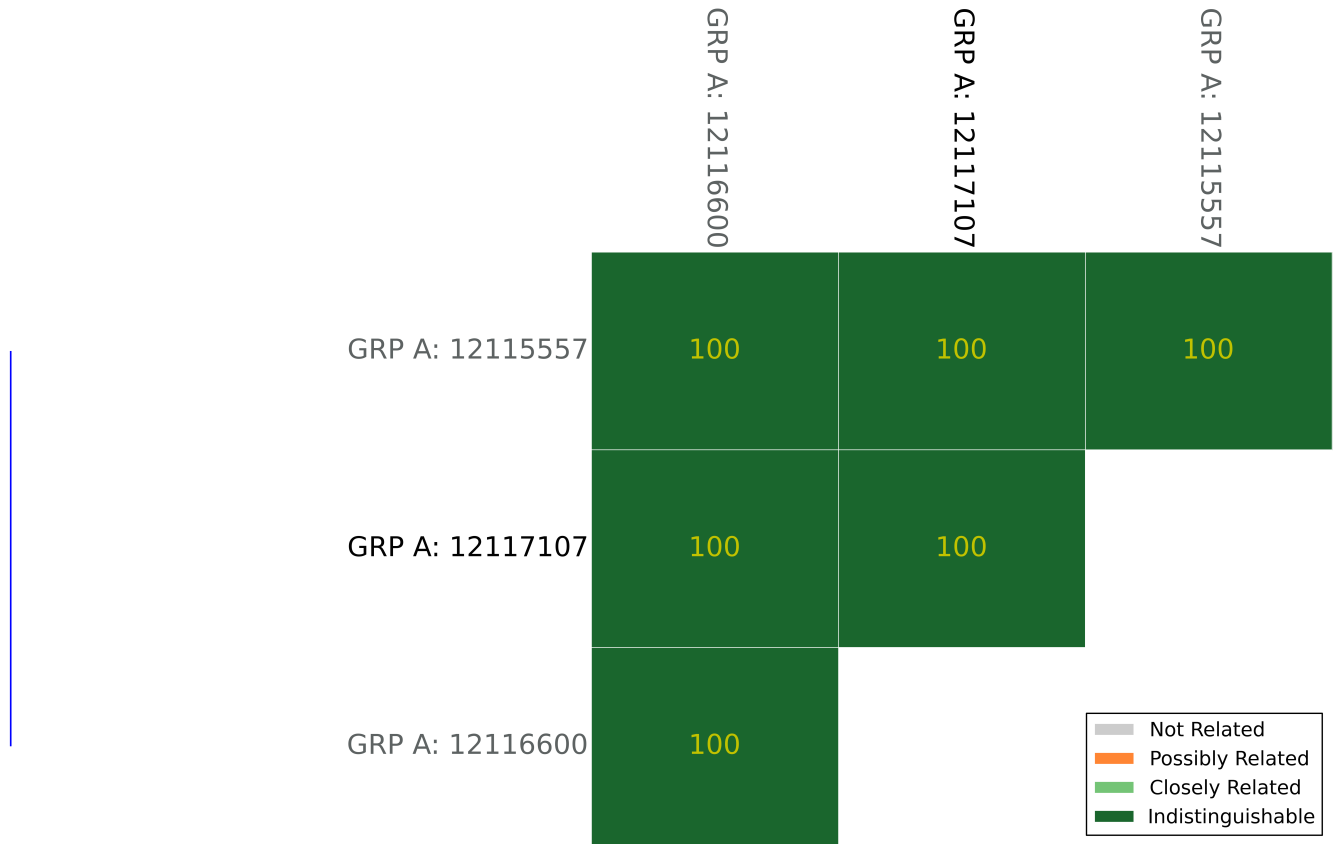

INTERPRETIVE DATA: Bacterial Strain Typing by NGS

Method

Whole Genome Sequencing (WGS) is performed using Ion Torrent sequencing chemistry. Reference-free pairwise comparisons are performed using short, overlapping sequence matching (kmer) analysis. Relationships are determined by the percent of kmers that match between isolate pairs.

Interpretation

Predicted relatedness is based on the total number of differences between the isolates, applying the thresholds shown in the table. The dendrogram and relationship matrix (see enhanced report) illustrate isolate relatedness. Interpretation of strain relatedness should be performed by an investigator knowledgeable about whole genome strain typing procedures and based on all available epidemiological evidence. Inferred relationships based on any strain typing method should not be used for individual patient management.

WGS Strain Typing provides substantial improvements in resolution and reproducibility when compared to pulsed-field gel electrophoresis (PFGE) and can be performed on a broad range of microorganisms. Test was validated for Staphylococcus, Acinetobacter, Enterococcus, Escherichia, Pseudomonas, Stenotrophomonas, Serratia, and Klebsiella species.




Patient:
ARUP Accession: 21-364-101995

Bacterial Strain Typing by Next Generation Sequencing

Patient: | Date of Birth: | Sex: | Physician:
Patient Identifiers: | Visit Number (FIN):

Category	Kmer Identity	Epidemiological Interpretation
Indistinguishable	≥ 99.9	Part of the outbreak
Closely related	99.8 - 99.2	Probably part of the outbreak
Possibly related	99.1 - 95.0	Possibly part of the outbreak
Unrelated	< 95	Not part of the outbreak

This test was developed and its performance characteristics determined by ARUP Laboratories. It has not been cleared or approved by the U.S. Food and Drug Administration. This test was performed in a CLIA certified laboratory and is intended for clinical purposes.

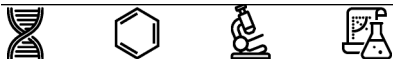

Patient:
ARUP Accession: 21-364-101995